



Thyroid Function	Fertility Hormones	Tumor Markers	Infectious Diseases*
TSH TT3 TT4 FT3 FT4 Anti-TG Anti-TPO TRAb* Tg*	β-HCG AMH E2 LH Prog PRL FSH Tes hGH 17-OH Prog* DHEA-SO4*	AFP CEA CA 125 CA 19-9 f-PSA t-PSA NSE* HE4* SCC* CYFRA21-1* proGRP* CA72-4* CA50* CA 15-3*	HBsAg HBsAb* HBeAg* HBeAb* HBcAb* HIV Ag/Ab anti-TP anti-HCV anti-HAV IgM* anti-HBc IgM* HBV PreS1* anti-HDV IgM* anti-HDV IgG* anti-HEV IgM* anti-HEV IgG*
Autoimmune ANA*	Torch Eugenics*	Hypertension*	Cardiac Markers
SSA(Ro52)* SSA(Ro60)* SSB(La)* Sm* Sm/RNP* Scl-70* Jo-1* His* AMA-M2* CENP-B* Rib-P* dsDNA* Nuc* PCNA* PM-Scl*	TOXO IgG* TOXO IgM* Rubella IgG* Rubella IgM* CMV IgG* CMV IgM* HSV 2 IgG* HSV 1/2 IgG* HSV 1/2 IgM*	ALD* Renin* Corti* ACTH* Angio-2*	MYO CKMB cTnI NT-proBNP
Prenatal Screening*	Autoimmune Diabetes*	Autoimmune Hepatitis*	
uE3* Free β-HCG* PAPP-A*	IAA* GADA* ICA* IA-2A*	AMA-M2* LKM* Sp100* SLA*	
Bone Metabolism	Glycometabolism	Inflammation	
25-OH Vit. D iPTH	C-Peptide Insulin	PCT IL6	
Anemia			
Ferritin Folate Vitamin B12			
Autoimmune Rheumatoid			
Anti-CCP			

The items marked with * are under validation

EN ISO 13485:2016




Toll Free No: 1800 425 7151

AGAPPE DIAGNOSTICS LTD.

“Agappe Hills”, Pattimattom (PO), Dist. Ernakulam, Kerala - 683 562, India.

TEL: + 91 484 2867000 | productcorp@agappe.in | www.agappe.com

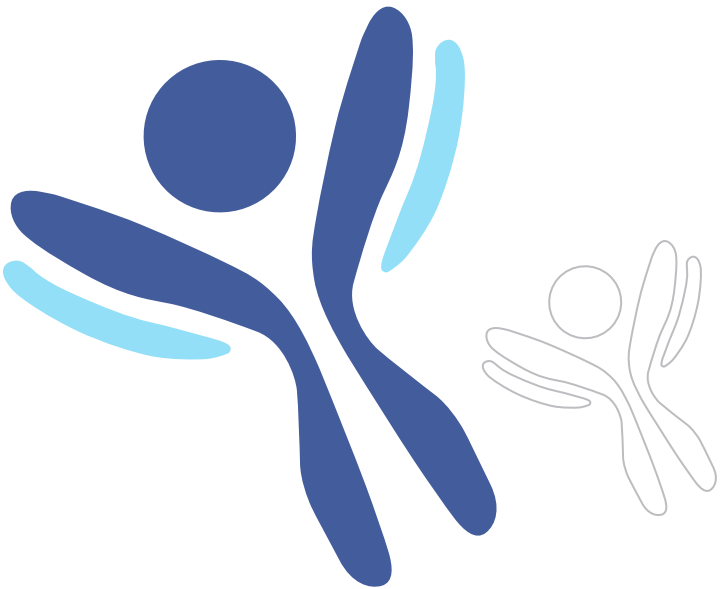

agappe

...an
IMMINENT
LOGIC

ADL/BB/Mispa-i200/R04-01/26 Company reserves the right to change any design and technical features of the product at any time, if needed.



Mispa i200

Automatic Chemiluminescence Analyzer

YOUR BEST PARTNER IN DIAGNOSTICS




Touch screen operation.
External mouse and key-
board can be connected



Used cuvettes recycling to
separate the solid and liquid,
which will efficiently reduce
the biological safety risks.




20 refrigerated
reagent positions,
temperature: 2-8 °C

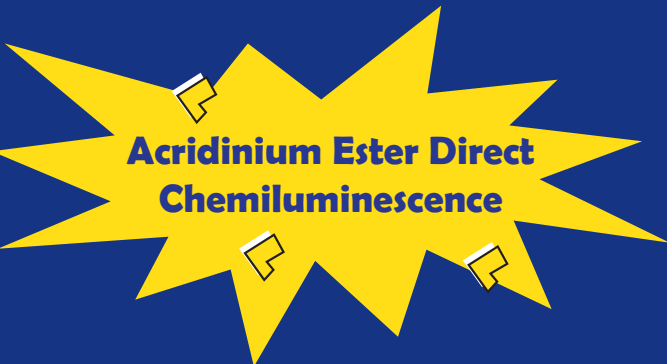


Ready to Use | Convenient Pack Sizes
CALIBRATORS INCLUDED IN THE KIT



Indian IVD Manufacturer to
Commission In-House CLIA
Reagents & Instrument
Production Line

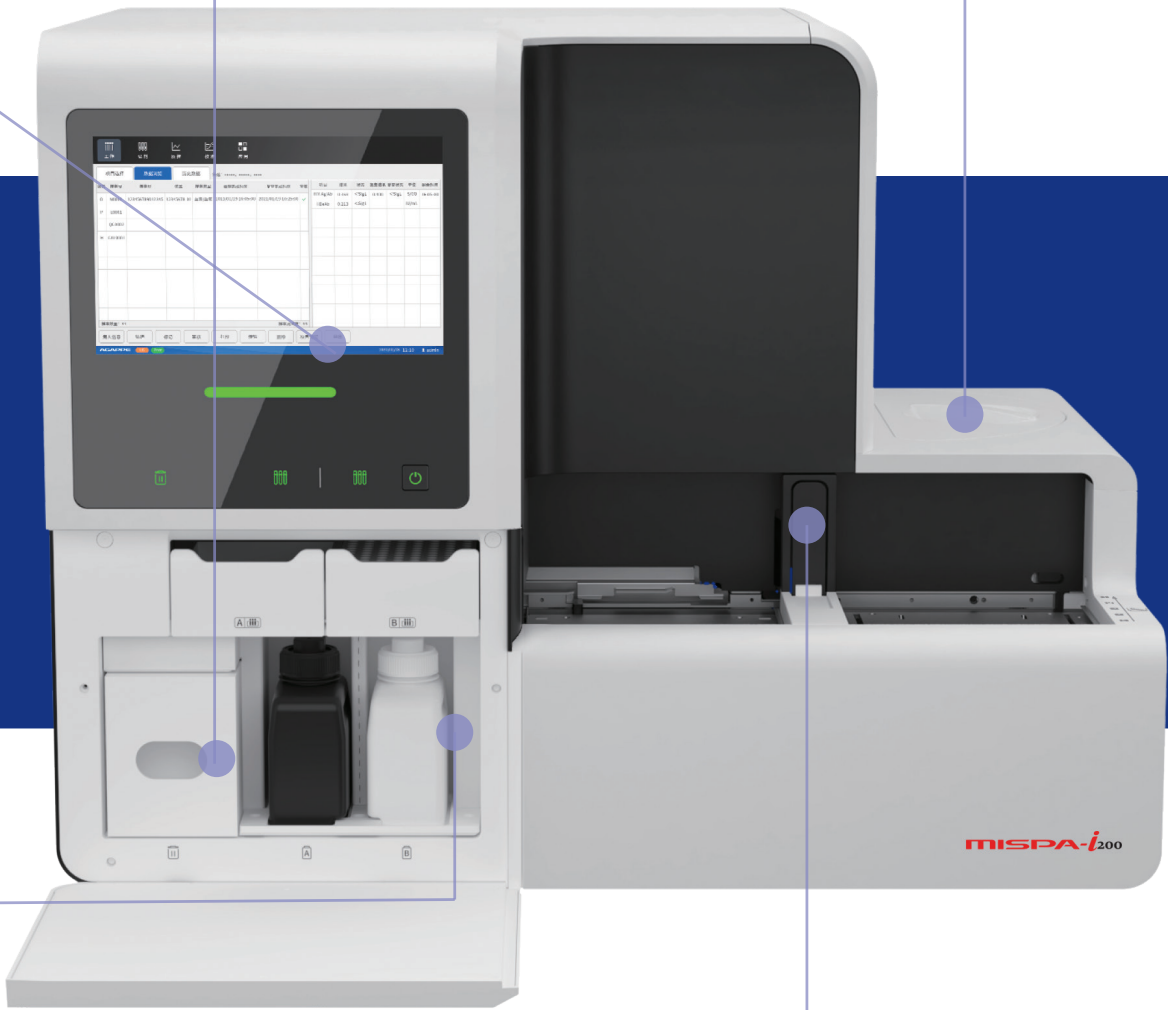

Activators with barcode,
special anti-reverse
design ensures correct
placement.

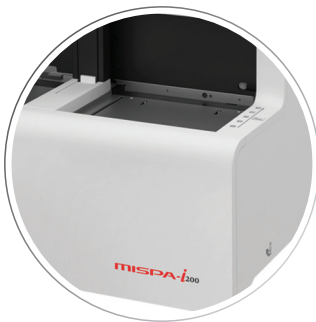

Dedicated STAT positions, STAT
sample can be detected emergen-

SPECIFICATIONS

Performance	Technology: Automatic /single tube Principle: Acridinium Ester (AE) Direct Chemiluminescence Throughput: 200 Tests/Hr Time to the first result: 14min for any cardiac profile test, 28 min for others
Sample Module	Sample loading capacity: 50, loading with sample racks. Automatic retest, STAT samples could be placed randomly and tested preferentially
Sample Loading System Module	Sample carry-over <0.1ppm Foam/Clot and Liquid level detection / Volume tracking / Collision-avoidance Sample type: Serum /Plasma /Whole blood /Urine
Reagent Module	20 reagent positions Information management through Radio Frequency Identification Card (RFID) Real-time, onboard mixture of magnetic beads Reagent online replacement
Reagent Refrigerator	Temperature: 2-8 °C
Detection Module	Reaction temperature: 37°C Mode: Compatible with the common immunoreaction Mixing: Non-contact Stir Magnetic separation: Multiple magnetic separation technology
Function	Operation system: Windows 7 and above Display: LCD touch screen Remote diagnosis system Automatic startup & shutdown Bi-directional communication

Environmental Requirements	Temperature: 10°C-30°C Humidity: 30%-70% Altitude: -200m-3700m
-------------------------------	--